

INVACARE®

Homecare Beds

- Beds
- Mattresses
- Bed Rails



Yes, you can.™

Invacare® Home Care Bed Series

Invacare® Home Care Beds deliver long-term savings over the lifecycle of the bed. The beds are durable and washable, reducing repair and cleaning costs, while their interchangeability reduces inventory. All of the home care bed features have been designed to make the bed easier to deliver, set up, clean and maintain.



Features

- **Bed components** – patented universal bed ends, headspring and footspring – are interchangeable with the original Invacare® Home Care Bed
- **Patented Universal bed ends simplify inventory** – store just one universal bed end instead of two different bed ends
- **Patented Universal bed ends** can be used as either the headboard or footboard on any Invacare bed (semi-electric, full-electric or manual, IVC or original) with the use of the IVC drive shaft with silver decal
- **High-impact bed end panels** are 22 times stronger than fiberboard and have better impact and scratch resistance
- **DC low voltage motor system** is faster, quieter and lighter than traditional AC systems
- **Motor system, bed frame, hand pendant and bed ends are completely washable** with standard water pressure (power washing is not recommended); motor system with waterproof enclosures is rated IP66
- **Ergonomic hand pendant comes with durable strain relief, heavy-duty cord** and an optional hand pendant cradle for convenient storage of hand pendant (part no. 26975)
- **Two locking and two non-locking casters** with non-marring rubber tread for secure positioning
- **Enhanced caster sockets** provide better retention during delivery and setup
- **Color-coded connections** on junction box make maintenance easier
- **Double-insulated junction box** eliminates need for grounding prong
- **Emergency crank** provided allows all motors to be operated manually in the event of a power failure
- **Extended maintenance schedule** – every three years or with every new user
- 350 lb. patient weight capacity
- 450 lb. total weight capacity
- UL962 and CSA certified

Invacare® Semi-Electric, Full-Electric and Manual Beds



Invacare® Semi-Electric Bed

Model no. 5310IVC

HCPCS code: E0295

The Invacare® Semi-Electric Bed combines effortless positioning of the upper body and knees with the economy of manual bed height adjustment. The pendant control provides motorized positioning of the upper body and/or knees. A manual crank located at the foot end of the bed adjusts bed frame height for caregiver assistance.



Invacare® Full-Electric Bed

Model no. 5410IVC

HCPCS code: E0297

The Invacare® Full-Electric Bed offers the greatest convenience for the patient and caregiver. The easy-to-use pendant control provides motorized positioning of the upper body and knees, and also adjusts bed frame height for efficient use of accessories and for safe and easy transfer.



Invacare® Manual Bed

Model no. 5307IVC

HCPCS code: E0293

The Invacare® Manual Bed offers the same convenient features as electric models, but adjusts manually for greatest economy. Manual cranks located at the foot end of the bed raise and lower upper body and knees, as well as bed height.

Specifications

5310IVC Semi-Electric Bed

Dimensions

Bed height*	15" min. - 23" max.
Sleep surface	80" L - 36" W
Overall bed	88" L - 36" W

	Carton size	Shipping weight
5000IVC Headspring	36" L x 44" W x 8" D	45 lb.
5890IVC Footspring	37" L x 47" W x 8" D	72 lb.
5301IVC Universal bed ends	25" L x 39" W x 7" D	55 lb.

Limited warranty

All welds	5 years
Mechanical/electrical	2 years

*Measurement from floor to bed deck.

Specifications

5410IVC Full-Electric Bed

Dimensions

Bed height*	15" min. - 23" max.
Sleep surface	80" L - 36" W
Overall bed	88" L - 36" W

	Carton size	Shipping weight
5000IVC Headspring	36" L x 44" W x 8" D	45 lb.
5490IVC Footspring	37" L x 47" W x 8" D	77 lb.
5301IVC Universal bed ends	25" L x 39" W x 7" D	55 lb.

Limited warranty

All welds	5 years
Mechanical/electrical	2 years

*Measurement from floor to bed deck.

Specifications

5307IVC Manual Bed

Dimensions

Bed height*	15" min. - 23" max.
Sleep surface	80" L - 36" W
Overall bed	88" L - 36" W

	Carton size	Shipping weight
5000IVC Headspring	36" L x 44" W x 8" D	45 lb.
5770IVC Footspring	37" L x 47" W x 7.5" D	62 lb.
5301IVC Universal bed ends	25" L x 39" W x 7" D	55 lb.

Limited warranty

All welds	5 years
Mechanical/electrical	2 years

*Measurement from floor to bed deck.

Whenever bed rails are used as part of any hospital or home care bed system, the potential for bed rail entrapment exists. To reduce this risk, a patient must always be assessed regarding their risk of entrapment before any product is used. Routine patient monitoring is also required once a product is in use. To learn more about the risk of bed rail entrapment please visit the US FDA's website at <http://www.fda.gov> and search for "Bed Rail Entrapment."

Invacare® Full-Electric Low and Value Care Semi-Electric Beds



Invacare® Full-Electric Low Bed

Model no. 5410LOW

HCPCS code: E0297

Invacare brings the safety and convenience of a low bed into the Home Care Bed Series with the Invacare® Full-Electric Low Bed. The split-spring design incorporated into the standard low bed profile combines the easy, one-person delivery of an Invacare® bed with the benefits of a low bed, where rails are not desirable or appropriate, but potential injuries from falling out of bed are still a concern. The Low Bed ends are fully compatible with existing Invacare full-electric bed frame inventories and the drive shaft, simplifying inventory management for providers.



Invacare® Value Care Semi-Electric Bed

Model no. VC5310

HCPCS code: E0295

The Invacare® Value Care Semi-Electric bed utilizes a single motor design making the bed lightweight and easy to assemble. It features the popular patented universal bed ends with the Invacare® deck design and a 9V battery backup for low cost protection from power outages. The advanced pendant design allows the head and knee sections to move either simultaneously or independently at the push of a button.

Specifications

Model no.	Bed deck height*	Overall bed	Sleep surface	Weight capacity	Limited warranty
VC5310	15" - 23"	88" L x 36" W	80" L x 36" W	450 lb.**	2 years

*Measured from floor to deck.

**Inclusive of all accessories.

Specifications

5410LOW Full-Electric Low Bed

Bed height*	9.5" min. - 20" max.	
Sleep surface	80" L - 36" W	
Overall bed	88" L - 36" W	
	Carton size	Shipping weight
5000IVC Headspring	36" L x 44" W x 8" D	45 lb.
5490IVC Footspring	37" L x 47" W x 8" D	77 lb.
5301LOW Universal bed ends	25" L x 39" W x 7" D	55 lb.

Limited warranty

All welds	5 years
Mechanical/electrical	2 years

*Measurement from floor to bed deck.



Features

Model no. VC5310

Bed Packages: VCPKGIVC-1633:

VC5310, 5185, 6629

VCPKGIVC2-1633:

VC5310, 5185, 6630

- Single motor design makes the bed lightweight and easy to assemble
- 9V battery backup provides convenient low cost protection in case of power outages
- Advanced pendant design allows head and knee sections to operate simultaneously or independently
- Patented universal bed ends are high impact and resist impacts and scratches
- Semi-electric operation:
 - Electric head and leg elevation
 - Manual hi/lo operation
- 5185 mattress (included in bed packages) meets 16CFR1633
- 350 lb. patient weight capacity
- 450 lb. total weight capacity



Invacare® BAR600 Bariatric Bed

The Invacare® Bariatric Bed is a heavy-duty full-electric bed frame designed for bariatric individuals. It is capable of supporting patients who weigh up to 600 lb.. The bed frame extends to the edges of the bed deck for better support than most other split-spring bariatric beds. Using the hand pendant, the consumer or caregiver is able to change the positioning of the head and foot sections, as well as the bed height. Despite being substantial in size, bed setup within the home is easy – no special tools are needed.



Hand pendant.



Invacare® Bariatric Bed

Model no. BAR600IVC

HCPCS code: E0301

The Invacare® BAR600IVC Bariatric Bed is a heavy-duty bed frame designed for bariatric individuals who weigh up to 600 lb. The easy-to-use pendant control provides motorized positioning of the head and foot sections, as well as the bed height.

Features

- Heavy-duty frame design ensures added strength for client support
- Heavy-duty bed sleep surface is 17% larger than conventional beds to ensure comfortable positioning
- Split-spring design with universal bed ends for ease in setup
- No special tools necessary for in-home setup
- Easy-release motor mounts allow quick motor removal
- Robotically welded, steel frame construction
- Pendant allows for dynamic client positioning
- Quiet, smooth operation
- Emergency crank manually functions to lower any raised section of bed
- Use with Invacare bariatric mattress (BARMATT42) and bed rails (BAR6640IVC)
 - BARMATT42 meets fire safety regulations 16CFR1633 and 16CFR1632
- 600 lb. client weight capacity

Specifications

BAR600IVC Bariatric Bed

Dimensions

Bed height*	min. 14" - max. 22"
Sleep surface length	80"
Sleep surface width	42"
Overall length	88"
Overall width	42"
Universal bed end	min. 27" - max. 34"

	Carton size	Shipping weight
BAR5000IVC Headspring	43.75" L x 10.5" W x 45" D	114 lb.
BAR5490IVC Footspring	47.25" L x 10.5" W x 45" D	97 lb.
BAR5301IVC Bed ends	45" L x 7.5" W x 23" D	64 lb.

Limited warranty

All welds	5 years
Frame/structure	5 years
Mechanical/electrical	1 year

*Measurement from floor to bed deck.

Invacare® BAR750 Bariatric Bed

The Invacare® BAR750 Bariatric Bed is a heavy-duty full-electric bed frame capable of supporting up to 750 lb. (650 lb. patient weight capacity). With the consumer in mind, this bed expands from 39" wide to a full 48" and from 80" in length to 88". Four ultra-quiet DC actuators easily reposition the head and foot sections, as well as the bed height. And, for the ultimate in safety, the bed comes standard with a battery backup.

HCPCS code: E0302



Features

- Heavy-duty frame design ensures added strength for client support
- Split-spring design eases delivery and setup
- Two locking casters on head section and two steering lock casters on foot section
- Two high-low motors, one on head section and one on foot section
- Pendant allows for dynamic client positioning
- Quiet, smooth operation
- Integrated battery backup, bed ends and half rails standard
- Choose from three mattress width options
- BARMATT Mattress options meet fire safety regulations 16CFR1633 and 16CFR1632
- 650 lb. patient weight capacity
- 750 lb. total weight capacity

Specifications

BAR750 Bariatric Bed

Dimensions

Bed height*	min. 15.5" - max. 29.5"
Sleep surface length	80" - 88"
Sleep surface width	39", 42" or 48"
Shipping dimension	56" L x 44" W x 41" D
Shipping weight	465 lb.

Limited warranty

Frame/structure	5 years
Mechanical/electrical	1 year

Model number

BAR750 Bariatric Bed, with 1 pair half rails

Bed Packages

BARPKG750-1633	BAR750/BARMATT39
BARPKG750-1-1633	BAR750/BARMATT42
BARPKG750-2-1633	BAR750/BARMATTEXP

Accessories

BARMATT39	Bariatric Mattress, 700 lb. capacity	80" L x 39" W x 6.5" D
BARMATT42	Bariatric Mattress, 700 lb. capacity	80" L x 42" W x 6.5" D
BARMATT48	Bariatric Mattress, 700 lb. capacity	80" L x 48" W x 6.5" D
BARMATTEXP	Bariatric Mattress, 700 lb. capacity	80" x 39" - 48" W x 6.5" D (expandable)
BAR750BOLSTER	Bed extension from 80" to 88"	8" x 39" (expandable)

*Measurement from floor to bed deck.

Invacare® Bed Rails

All Invacare® Bed Rails feature durable chrome-plated, welded steel construction and are designed to withstand heavy-duty use. These rails can be installed permanently or semi-permanently on most spring-fabric medical-style beds. When properly used, bed rails are some of the most important safety and accident prevention accessories used in the patient room today.



Invacare® Reduced Gap Full-Length Bed Rail

HCPCS code: E0310

Model no. 6629

- Durable 7/8" chrome-plated, welded steel construction
- Self-adjusting spring-loaded crossbraces snap easily into place without tools for quick installation



Invacare® Reduced Gap Half-Length Bed Rail

HCPCS code: E0305

Model no. 6630

- Heavy-duty 1" chrome-plated, welded steel construction
- Dual flat-channel crossbraces are secured with wing-nut clamp assemblies (some assembly required)
- For optimal safety and convenience, four half-length rails are recommended (one set for head area and one set for foot area)



Invacare® Assist Bed Rail

HCPCS code: E0305

Model no. 6632

- Heavy-duty 1" powder coated, welded steel construction
- 3 positions: Guard, Assist, Transfers
- Moves easily out of the way with a quick release button



Invacare® Extra-Wide Full-Length Bed Rail

HCPCS code: E0310

Model no. 6628-TSS

- Heavy-duty 1" chrome-plated, welded steel construction
- Self-adjusting extra-wide spring-loaded crossbraces snap easily into place without tools for quick installation
- Optional hardware is included for clamp-on
- For use with TSS mattresses



Invacare® Reduced Gap Heavy-Duty Half-Length Bed Rail

HCPCS code: E0305

Model no. BAR6640IVC

- Heavy-duty 1" chrome-plated, welded steel construction
- Bolt-on hardware for added strength
- For optimal safety and convenience, four half-length rails are recommended (one set for head area and one set for foot area)
- Use with BAR600IVC bed only

Specifications

Model no.	Rail length*	Rail height above deck		Width between bed rails	Shipping weight	Limited warranty
		Raised	Lowered			
6628-TSS	56"	14.5"	3.75"	38"	31.5 lb.	1 year
6629	55"	14.5"	4.25"	36"	27.5 lb.	1 year
6630	32"	15"	4.75"	36.5"	23.5 lb.	1 year
6632	27"	14.5"	n/a	36"	21 lb.	1 year
BAR6640IVC	28.75"	15"	3.50"	42"	25 lb.	1 year

*Measurement on flat bed with crossbraces installed as per instruction sheet.

Whenever bed rails are used as part of any hospital or home care bed system, the potential for bed rail entrapment exists. To reduce this risk, a patient must always be assessed regarding their risk of entrapment before any product is used. Routine patient monitoring is also required once a product is in use. To learn more about the risk of bed rail entrapment please visit the US FDA's website at <http://www.fda.gov> and search for "Bed Rail Entrapment."

Invacare® Mattresses, Mattress Covers

Now you can offer your customers the comfort, protection and durability of superior home care mattresses by Invacare, along with Invacare® mattress covers. Choose from innerspring or foam mattresses. Mattress covers will keep your Invacare® mattress clean and fresh, and protect it from stains, odors, dust and incontinence. Covers are made of non-allergenic waterproof vinyl and are available in zippered and contoured styles.



Invacare® Innerspring Mattress

Model no. 5185, 5185XL

HCPCS code: E0271

- High quality .312 innerspring design
- Durable 11-ounce vinyl cover is anti-bacterial, anti-static, acid-resistant and waterproof for easy use and care
- 5185XL (not shown) is available for use with extender kit
- Meets fire safety regulations 16CFR1633 and 16CFR1632



Invacare® MC0195-1 zippered mattress cover comes individually wrapped; 12 per case.

Invacare® Economy Foam Mattress

Model no. 5180, 5184

HCPCS code: E0272

- High quality 1.5 lb. density foam provides comfort and durability in a one-piece construction
- Meets fire safety regulations 16CFR1633 and 16CFR1632



Invacare® MC0065-1 contoured mattress cover comes individually wrapped; 12 per case.



Invacare® Bariatric Foam Mattress

Model no. BARMATT39, BARMATT42, BARMATT48, BARMATTEXP

HCPCS code: E0272

- Constructed of latex-free materials
- Durable 11-ounce vinyl cover is anti-bacterial and waterproof for easy use and care
- Solid inner core
- Mattress is reversible and flippable
- Constructed of high resilience foam for extended life and service for the homecare provider
- Easy to clean
- Meets fire safety regulations 16CFR1633 and 16CFR1632

Specifications

Model no.	MC0165-1	MC0095-1	5180	5184	5185	5185XL	BARMATT39	BARMATT42	BARMATT48	BARMATTEXP
Dimensions										
Overall length	80"	80"	80"	84"	80"	84"	80"	80"	80"	80"
Overall depth	6"	6"	5.5"	5.5"	6"	6"	6.5"	6.5"	6.5"	6.5"
Overall width	36"	36"	36"	36"	36"	36"	39"	42"	48"	39"-48" expandable
Gauge	4"	4"	n/a	n/a	n/a	n/a	n/a	n/a	n/a	n/a
Limited warranty	30 days	30 days	1 year	1 year	1 year	1 year	1 year	1 year	1 year	1 year

Invacare® Group I Therapeutic Support Surfaces

Invacare® Group I Therapeutic Support Surfaces are designed to aid in the prevention and treatment of pressure ulcers. Product types include an alternating pressure pad system, gel overlay and a foam mattress replacement. All of these products are simple to use and easy to maintain.



Invacare® Solace® Prevention 0080 Foam Mattress

Model no. SPS0080

HCPCS code: E0184

- Zoned mattress relieves pressure in the important susceptible body zones
- Helps distribute weight and lower interface pressure to help reduce the risk of pressure ulcers in low risk patients
- Vinyl non-slip bottom cover provides secure positioning on bed sleep surfaces
- Vapor permeable, fluid-resistant cover with a 70 denier nylon top to reduce friction
- Die cut foam helps circulate air around mattress user and facilitates gatching
- Single layer quality foam designed for reducing the risk of pressure ulcers
- Reversible head-to-foot for even wear
- Available in 80" and 84" lengths
- Does not contain latex
- Meets fire safety regulations 16CFR1633 and 16CFR1632



Invacare® Gel Foam Mattress Overlay

Model no. IVCGFM02

HCPCS code: E0185

- Reversible for increased service life
- Durable aggressive base material with strong elastic corner securing straps reduces overlay movement
- Vapor permeable, fluid-resistant cover with low friction and shear (do not machine wash)
- Meets California Technical Bulletin #117 for fire retardancy
- Does not contain latex
- .36mm gauge of vinyl



Invacare® CareGuard™ APP Alternating Pressure Pad System

Model no. CG9701

HCPCS code: E0181

- Advanced pump technology alternately inflates and deflates air cells
- Variable pressure setting for added patient comfort
- Pad includes flaps that wrap around mattress for better stability
- Pump can be mounted on the bed with built-in bracket
- Pad meets California Technical Bulletin #116 for fire retardancy
- Uses CG9702 CareGuard pads only
- Does not contain latex

Specifications

Model no.	Weight	Length	Width	Thickness	Power	Code approval	Weight capacity	Limited warranty
SPS0080	19 lb.	80"	36"	6"	n/a	n/a	350 lb.	2 years (foam only) 1 year on cover
IVCGFM02	55 lb.	79"	35"	3.5"	n/a	n/a	250 lb.	6 months
CG9701	Pump: 3.5 lb. Pad: 4.6 lb.	10" 118"	6" 34"	4" 2.5"	115V, 160 HZ	UL-1431	250 lb.	1 year 30 days
CG9702	4.6 lb.	118"	34"	2.5"	n/a	n/a	250 lb.	30 days

Invacare® Trapeze Equipment


The Invacare® trapeze unit is an important patient room accessory designed to help the users change positions while in bed and to aid in transferring from bed to chair with easy attendant assistance. Trapeze units have a wide range of height adjustments and handbar positions to maximize patient accessibility.

Invacare® Trapeze Floor Stand

Model no. 7714P

HCPCS code: E0940 (only for base)

- Provides free-standing support for model no. 7740P fixed trapeze unit
- Folds flat for convenient storage and handling
- Assisted weight capacity of 250 lb.

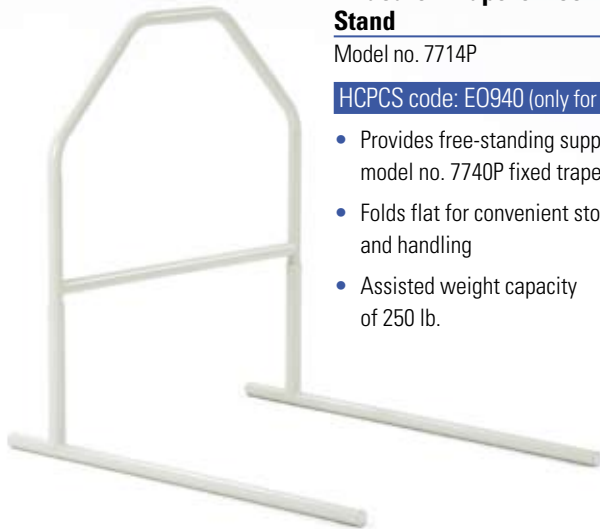

Invacare® Fixed Offset Trapeze Bar with Trapeze

Model no. 7740P

HCPCS code: E0910

- Heavy-gauge, octagon steel construction; no tools required for easy assembly
- Wall bumpers and padded mounting brackets protect bed ends and wall surfaces
- For use on model no. 7714P floor stand or on any Invacare home care bed head section
- Assisted weight capacity of 250 lb.
- Includes bar, brackets and grab bar



7714P Shown with 7740P Trapeze

Specifications

Model no.	7714P	7740P
Dimensions		
Overall height	38"	61"
Octagon steel tubing, diameter	1.25"	1.25"
Reach	n/a	34"
Base width	32"	n/a
Base length	40"	n/a
Product weight	24 lb.	17 lb.
Shipping	via UPS®	via truck
Limited warranty	3 years	3 years

Note: Do not attach trapezes to the end of the bed. Trapezes may be attached to the head section of the bed end only.



Invacare® Bariatric Floor Stand with Trapeze

Model no. BARTRAP

HCPCS code: E0912

- Three locking overhead boom angle positions for comfortable, secure use
- Tool-free design simplifies assembly
- Four piece sub-assembly significantly reduces setup time
- Lightweight design for easy delivery
- Free-standing structure permits use with Bariatric, Home Care or non-medical bed
- Assisted 1000 lb. patient weight capacity



Specifications

Model no.	BARTRAP
Dimensions	
Overall height	73"
Reach	35.5"
Base width	40"
Base length	50.5"
Product weight	72 lb.
Shipping weight	80 lb.
Limited warranty	3 years

Invacare® Bed Accessories

Invacare® overbed tables are practical, easy-to-use units suitable for both home and institution. The table height easily adjusts for comfort while reading, writing or eating. These tables feature an attractive walnut wood-grained laminate top with a T-molded edge that is designed to resist damage. The heavy gauge tubular welded steel construction ensures durability, and the 1.5" casters ensure smooth mobility of the table in any direction.



Invacare® Auto-Touch Overbed Table

Model no. 6417

- Light upward touch on patented spring-loaded locking handle allows table to elevate freely
- Light downward touch on table top locks table in position

Invacare® Tilt-Top Overbed Table

Model no. 6418

- Adaptable tilt-top allows three angle positions in either direction using tilt-release lever
- Light upward touch allows table to elevate freely
- Light downward touch on table top locks it into position

Specifications

Model no.	6418	6417
Dimensions		
Table top	30" L x 15" W x .75" D	30" L x 15" W x .75" D
Base	26" L x 15.5" W	26" L x 15.5" W
Height adjustment	25.75" - 39"	29" - 45"
Weight	28 lb.	27 lb.
Column/base/ support bracket	14/16/13 gauge	14/16/13 gauge
Weight capacity	25lb.	25lb.
Limited warranty	1 year	1 year



Bed Buggy

Model no. 5140

- Used to transport split-spring or assembled bedspring quickly and easily
- Allows for one-person setup in much less time
- Can be used on all Invacare beds



Extender Kit

Model no. 5143IVC for Invacare Home Care Beds

- Designed to expand 80" bed frame to 84"
- Complete with extended head pull tube and drive shaft, link fabric and helicals
- Not available for bariatric bed
Model no. BAR600IVC



Stor-A-Bed

Model no. 5142

- The ultimate bed storage unit. Uses both bed ends to transport or store bed. Helps maintain bed inventory, requiring less warehouse space
- Can be used on all Invacare beds



Extender for Value Care Beds

Model no. 5143A

- Designed to expand 80" bed frame to 84" bed frame
- Attaches at the head end of the bed
- Works on the ValueCare bed (Item Number: VC5310)
- Complete with frame extension, sleep surface extension, and extended drive shaft

HCPC Guide for Invacare® Home Care Beds and Bed Packages

Full-Electric Beds and Bed Packages

Model Number	Description	HCPCS Code	HCPCS Description
5410IVC 5410LOW	5410IVC 5410LOW	E0297	Hospital bed, total electric (head, foot and height adjustments), without side rails, without mattress.
BED2-1633 BED9-1633 BED33-1633 BED37-1633 BED38LOW-1633 BED39LOW-1633 BED21-1633 BED23-1633 BED44LOW	5410IVC, 6629, 5185 5410IVC, 6630, 5185 5410IVC, 6629, 5180 5410IVC, 6630, 5180 5410LOW, 6632, 5185 5410LOW, 6632, 5180 5410IVC, 6629, SPS0080 5410IVC, 6630, SPS0080 5410LOW, 6632, SPS0080	E0265	Hospital bed, total electric (head, foot and height adjustments), with any type side rails, with mattress.

Semi-Electric Beds and Bed Packages

Model Number	Description	HCPCS Code	HCPCS Description
5310IVC VC5310	5310IVC VC5310	E0295	Hospital bed, semi-electric (head and foot adjustment), without side rails, without mattress.
BED4-1633 BED10-1633 BED32-1633 BED36-1633 BED20-1633 BED22-1633 VCPKGIVC-1633 VCPKGIVC2-1633	5310IVC, 6629, 5185 5310IVC, 6630, 5185 5310IVC, 6629, 5180 5310IVC, 6630, 5180 5310IVC, 6629, SPS0080 5310IVC, 6630, SPS0080 VC5310, 6629, 5185 VC5310, 6630, 5185	E0260	Hospital bed, semi-electric (head and foot adjustment), with any type side rails, with mattress.

Bariatric Beds and Bed Packages

Model Number	Description	HCPCS Code	HCPCS Description
BAR600IVC	BAR600IVC	E0301	Hospital bed, heavy-duty, extra wide, with weight capacity greater than 350 lb., but less than or equal to 600 lb., with any type side rails, without mattress.
BAR750	BAR750 (Rails Standard)	E0302	Hospital bed, extra heavy-duty, extra-wide, with weight capacity greater than 600 lb., with any type side rails, without mattress.
BARPKGIVC-1633	BAR600IVC, BAR6640IVC, BARMATT42	E0303	Hospital bed, heavy-duty, extra-wide, with weight capacity greater than 350 lb., but less than or equal to 600 lb., with any type side rails, with mattress.
BARPKG750-1633 BARPKG750-1-1633 BARPKG750-2-1633	BAR750, BARMATT39 (Rails Standard) BAR750, BARMATT42 (Rails Standard) BAR750, BARMATTEXP (Rails Standard)	E0304	Hospital bed, extra heavy-duty, extra-wide, with weight capacity greater than 600 lb., with any type side rails, with mattress.

Manual Bed and Bed Packages

Model Number	Description	HCPCS Code	HCPCS Description
5307IVC	5307IVC	E0293	Hospital bed, variable height, hi-lo, without side rails, without mattress.
BED8-1633 BED12-1633	5307IVC, 6629, 5185 5307IVC, 6630, 5185	E0255	Hospital bed, variable height, hi-lo, with any type side rails, with mattress.

Invacare Corporation
www.invacare.com

USA

One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada

570 Matheson Blvd. E.
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

© 2011 Invacare Corporation.
All rights reserved.
Trademarks are identified by the symbols ™ and ®.
All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted.
UPS is a registered trademark of United Parcel Services of America, Inc. Specifications are subject to change without notice.

Form No. 00-276
Rev. 04/11
100803

All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further details.

