

INVACARE®  
**FULL-ELECTRIC LOW  
BED**

---



---

Invacare brings the safety and convenience of a low bed into the Homecare Bed Series with the **Invacare® Full-Electric Low Bed**. The split-spring design incorporated into the standard low bed profile combines the easy, one-person delivery of an Invacare bed with the benefits of a low bed, most often used when injuries from falling out of bed are a concern. The Low Bed ends are fully compatible with existing Invacare full-electric bed frame inventories and the drive shaft, simplifying inventory management for providers.

Invacare beds, mattresses and rails, excluding therapeutic support surfaces and bariatric beds, when used as a system and sold new and unused after August 1, 2007, meet the voluntary FDA Guidelines for bed rail entrapment. Visit [www.invacare.com/bedsafety](http://www.invacare.com/bedsafety) for additional safety and entrapment information.

---

## Features

High impact bed end panels are durable, impact and scratch resistant; washable for easy cleaning

Universal low bed end can be used as either the headboard or footboard on the Low Bed

Interchangeability simplifies storage, reduces inventory complexity

Rubber caster tread for better tracking with non-marring material

DC low voltage, compact motor system has quiet motor operation and fast cycle time on lowering height of bed for efficient bed set-up

Washable motor system and frame—Motor system with waterproof enclosures is rated IPX6. Washable with standard water pressure. The bed must not be power washed


Color-coded connections on junction box for easier maintenance

Double-insulated junction box helps provide extra protection against electrical shock

Ergonomic hand pendant comes with durable strain relief and heavy-duty cord. Waterproof casing for easy cleaning, conveniently attaches to bed sheet with clip or bed rail with hand pendant cradle (sold separately)

Hand pendant cradle available for convenient storage of hand pendant (model # 26975)

## HCPCS Code: E0297

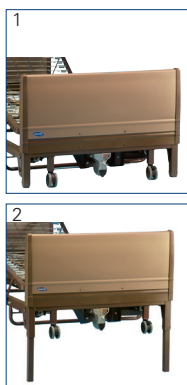

1. Bed in lowest position  
2. Bed in highest position



**Invacare Corporation**  
www.invacare.com

### USA

One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 333-6900

### Canada

570 Matheson Blvd. E.,  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

©2021

Invacare Corporation.  
All rights reserved.  
Trademarks are identified by the symbols ™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.  
Form No. 04-173  
Rev. 03/21  
210353



Model #	Dimensions	Sleep Surface	Weight Capacity	Shipping Weight	Carton Dimensions	Regulation Compliance	Limited Warranty
5410LOW	88" L 9.5"-20" H 36" W	80" L 36" W	Safe Working Load: 450 lb.* Max. User Weight: 350 lb.	N/A	N/A	UL962 CSA Certified	5 Years Welds 2 Years Mechanical/ Electrical
5000IVC Headspring				45 lb.	43" L 8" W 36" D		5 Years Welds 2 Years Mechanical/ Electrical
5490IVC Footspring				81 lb.	47" L 9" W 36" D		5 Years Welds 2 Years Mechanical/ Electrical
5301LOW Universal Bed End				55 lb.	39" L 7" W 24" D		

\*Includes user, mattress and all accessories.

All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.