


INVACARE®
TDX® SI
POWER WHEELCHAIR

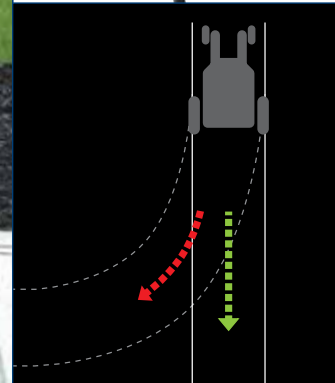

The **Invacare® TDX® SI Power Wheelchair** offers a streamlined approach to fully featured rehab power wheelchairs. With standard features such as Invacare® MK6i™ Electronics, simplified Invacare® SureStep® Suspension with Stability Assist and a top speed of 5 mph, the TDX SI Wheelchair delivers performance in an elegant design. The TDX SI Wheelchair offers several seating options including an adjustable rehab seat, captain's seating and Invacare® Formula™_{CG} Tilt. With a focus on excellent standard features, the TDX SI Wheelchair is the only name you need to know.


Features

Invacare® Formula™ CG Powered Tilt available
 Adjustable Rehab (ASBA) seat
 Semi-recline captain's seat available in 16", 18", 20" and 22" widths and 16"-20" depths
 Center-wheel drive and narrow base combine for tight turning radius and easy access through doors and hallways
 Integrated rear suspension provides a non-jarring, smooth ride when tackling rough terrain or thresholds
 Removable battery access door with slide-out battery tray makes servicing a snap
 Transport Ready Option (TRRO) available*

Invacare Core Technologies

Invacare® SureStep® Suspension with Stability Lock provides a smooth and stable transition
 Invacare® G-Trac™ Technology provides precise tracking for imprecise environments
 Invacare® MK6™ Electronics are the benchmark in innovative technology

HCPCS Code: K0848, K0849, K0856

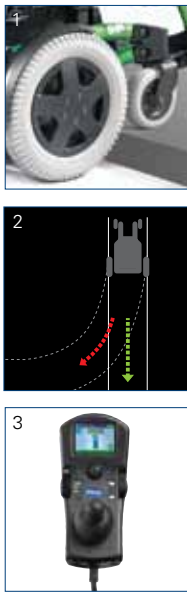

1. Simplified Invacare® SureStep® Suspension with Stability Assist
2. Invacare® G-Trac™ Technology
3. Invacare® MK6™ Electronics



*TRRO includes four factory-installed brackets and wheelchair anchored pelvic belt. TRRO has been crash-tested in accordance with ANSI/RESNA WC Vol 1. As of this date, the Department of Transportation has not approved any tie-down systems for transportation of a user while in a wheelchair, in a moving vehicle of any type. It is Invacare's position that users of wheelchairs should be transferred into appropriate seating in vehicles for transportation and use be made of the restraints available by the auto industry. Invacare cannot and does not recommend any wheelchair transportation systems.

Invacare Corporation
 www.invacare.com

USA
 One Invacare Way
 Elyria, Ohio
 44035-4190
 (800) 333-6900

Canada
 570 Matheson Blvd. E.,
 Unit 8
 Mississauga, Ontario
 L4Z 4G4 Canada
 (800) 668-5324

©2012
 Invacare Corporation.
 All rights reserved.
 Trademarks are identified by the symbols ™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.
 Form No. 07-073
 Rev. 09/12
 120336

Model #	Overall Dimensions**	Seat Dimensions	Seat-to-Floor Height	Drive Wheels	Ground Clearance	Turning Radius	Incline
TDXSI, TDXSI-CG, TDXSIV	40.6" L 25.4" W 35.5" H	12"-22"x 12"-20"	16.5"-21"	14"	3"	20.5"	9°
	Total Weight	Weight Capacity	Motor	Max Speed	Battery	Battery Range	Limited Warranty
	Up to 295 lb.	300 lb.	2-Pole	5 mph	22NF	Up To 17 Miles	5 Years Base Frame 3 Years Seat Frame 1 Year Electronics/Motor/Gearbox

**Rehab seat and 90° swingaway footrest
 Do not operate on roads, streets or highways.

All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.

