

INVACARE®
microAIR® MA80
SERIES THERAPEUTIC SUPPORT SURFACE



The **Invacare® microAIR® MA80 Therapeutic Support Surface** model is a true low air loss mattress that is ideal for situations where maceration and moisture management are an issue. The MA80 mattress helps manage skin moisture to avoid maceration and maintain a dryer environment for pressure ulcer healing in the treatment of Stage I to IV pressure ulcers.

The **Invacare® microAIR® MA85 Therapeutic Support Surface** model builds on the true low air loss model benefits, while increasing the therapy with the addition of alternating pressure to reduce surface interface and offer pressure relief.

Features

Blower – Provides greater air output, quick set up and more efficient therapy cycles

Wireless auto fowler – When the transmitter in the head section is raised, 25% more air is added into the center section of the bed, preventing the patient from bottoming out

Audio and visual alarms – Warn patient and caregiver of low pressure and power outage

Lock out button – Prevents against tampering or accidental adjustments
 Power outage safety – Returns to pre-power outage settings

Cover – Dual layer quilted cover, anti-microbial, fluid-proof, vapor permeable, guards against friction and denier nylon bottom shear


Settings – Based on patient weight to ensure correct pressure

Constant pressure monitoring – Efficient real time monitoring for correct pressure relief

Five pressure zones – Maintained automatically for softer head and heel sections and firmer torso pressure interfaces

Pulse setting – Available on MA80

HCPCS Code: E0277

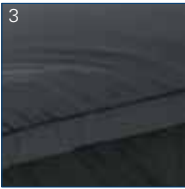

- 1. Interchangeable air cells
- 2. Audio and visual alarms
- 3. Stretch cover, denier nylon bottom



Invacare Corporation
 www.invacare.com

USA
 One Invacare Way
 Elyria, Ohio
 44035-4190
 (800) 333-6900

Canada
 570 Matheson Blvd. E.,
 Unit 8
 Mississauga, Ontario
 L4Z 4G4 Canada
 (800) 668-5324

©2012
 Invacare Corporation.
 All rights reserved.
 Trademarks are identified by the symbols ™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.
 Form No. 10-130
 Rev. 03/12
 111234



Model #	Dimensions	Weight Capacity	Blower	Alarms	Regulation Compliance	Limited Warranty
MA80	80" L 36" W 10" H	350 lb.	1275 LPM	Audio and Visual	16 CFR 1632 16 CFR 1633	1 Year (Cover/Mattress) 2 Years (Power Unit)
MA80B42/ MA8048	80" L 42"/48" W 10" H	1,000 lb.	1275 LPM	Audio and Visual	16 CFR 1632 16 CFR 1633	1 Year (Cover/Mattress) 2 Years (Power Unit)
MA85, MA85RSR Raised side rail option	80" L 36" W 10" H	350 lb.	1275 LPM	Audio and Visual	16 CFR 1632 16 CFR 1633	1 Year (Cover/Mattress) 2 Years (Power Unit)
MA85B42/ MA85B48	80" L 42"/48" W 10" H	1,000 lb.	1275 LPM	Audio and Visual	16 CFR 1632 16 CFR 1633	1 Year (Cover/Mattress) 2 Years (Power Unit)

All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.