

INVACARE®
RELIANT™ 600
HEAVY-DUTY LIFT



Designed to meet the special needs of your patients weighing up to 600 lb., the **Invacare® Reliant™ 600 Heavy-Duty Lift** is secure and sturdy, providing a safe alternative to manual lifting. Both durable and functional, the Reliant 600 Lift is also a reassuring and comfortable patient lift. When answering the needs of every patient is your priority, this is your wisest choice.

Invacare® Reliant™ 600 Lift Features

High lift range allows transfer to high surfaces
Low friction casters significantly improve rollability and contribute to caregiver safety
Designed to meet the needs of the growing bariatric population up to 600 lb.
Manual lowering for instant response to power loss
Unique sling fabric designed to provide the utmost in comfort and strength

Available with power opening base (model # RPL600-2)
Essential device for safe patient handling in hospitals and LTC (long-term care) facilities
Limited warranty includes three years on lift, one year on electric components
Rear locking casters are standard
Manufactured in an ISO13485 facility

HCPCS Code: E0635



RPL600-1: Powered with low base



RPL600-2: Powered with power opening low base



1. **Manual Lowering**

Manual emergency lowering gently and safely lowers the patient in case of total power loss

Essential safety feature

Manual lowering is easy and effortless. Simply insert the end of a ball point pen into the hole marked "emergency" or lift on the accessible ring-pull



2. **Invacare® Reliant™ Digital Scale**

Retrofits on all Reliant 450/600 lifts. Model # 9805 requires Reliant swivel bar

Easy to operate, four-button display

Automatic shut-off feature prolongs battery life

Two-year limited warranty



3. **Electronics**

Removable battery box conveniently located on mast features an emergency off button and beeps when battery needs recharged. Requires one 24-volt battery

Quick-release battery eliminates downtime

Anti-entrapment feature stops motion if boom meets resistance while lowering



4. **Charging Station**


Battery can be replaced immediately with second battery in wall-mounted charger

UL listed and CSA certified



1. **Mesh sling model # R140**

2. **Mesh sling with commode opening model # R141**

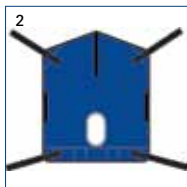

Invacare Corporation
www.invacare.com

USA
One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada
570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

©2012
Invacare Corporation.
All rights reserved.
Trademarks are identified
by the symbols ™ and
©. All trademarks are
owned by or licensed
to Invacare Corporation
unless otherwise noted.
Specifications are subject
to change without
notification.
Form No. 11-245
111281



Model #	Height At Sling Hook Up	Base Dimensions	Casters	Weight	Weight Capacity	Battery Life Indicator	Emergency Controls
RPL600-1	28"-68" H	4.5"H Clearance 41"W Open 26.5"W Closed 48"L	3.5" Front 5" Rear	112 lb.	600 lb.	N/A	Stop and Down
RPL600-2	28"-68" H	4.5"H Clearance 41"W Open 26.5"W Closed 48"L	3.5" Front 5" Rear	115 lb.	600 lb.	N/A	Stop and Down



	Battery	Charge Time	Charger Input	Charger Output	Limited Warranty
RPL600-1	N/A	N/A	N/A	N/A	3 Years Lift 1 Years Electronics
RPL600-2	N/A	N/A	N/A	N/A	3 Years Lift 1 Years Electronics