
INVACARE® |-L|FT™ SERIES LIFTS



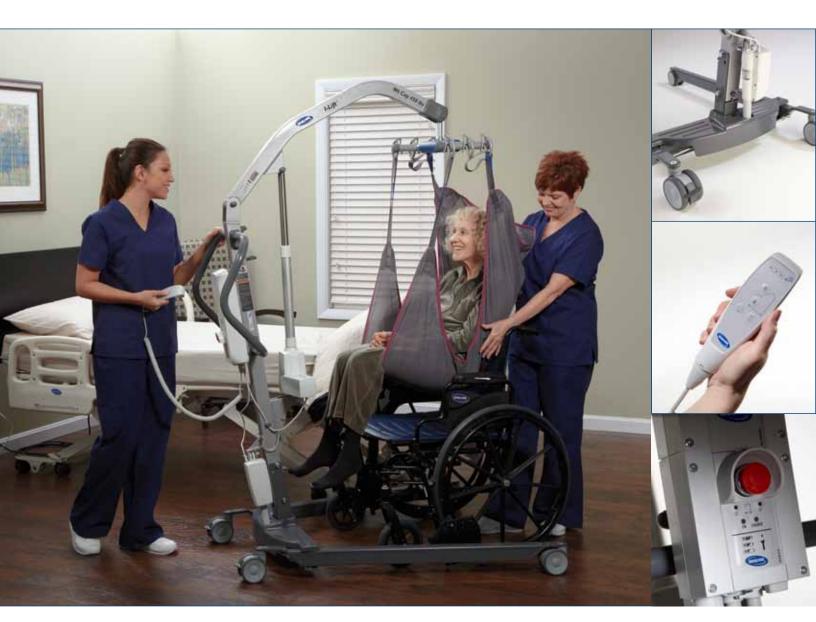

Raising the Bar for Safe Patient Handling Equipment

Designed to handle the most common transfer situations with ease, the Invacare® I-Lift™ Series with a 450 lb. capacity will accommodate patient transfers from a bed, a stretcher, a wheelchair, a commode and even from the floor.

With increased maneuverability, intuitive ease of use and enhanced safety features, the I-Lift Series lifts deliver advanced features that help protect and assist caregivers in providing safe and effective healthcare to their patients.

The Invacare® I-Lift™ Series was designed from the ground up to deliver the *safety, maneuverability* and *ease of use* you expect from a premium lift.

Safety

Resident and Caregiver safety is in the details. I-Lift Series Lifts were designed with innovative features to meet or exceed existing international hoist standards (ISOEN10535).

Sleek European styling with no exposed pinch points

Integrated anti-entrapment features automatically stops the lift if it meets any kind of resistance while it is being lowered

Large aluminum safety clips provide confidence that slings are firmly attached

Maneuverability

Caregivers will appreciate unsurpassed maneuverability in tight spaces with a low base clearance and small turning radius

Ergonomic dual angled push handle better accommodates a variety of caregivers' heights and allows access from both sides

Premium casters and base design require 60% LESS force to push/pull the lift

Ease of Use

Wide base opening is effortless with the push of a button

Powder coated, protective leg and base cover make the lift easy to clean and disinfect

Both models include

- A slim, lightweight and ergonomic hand pendant with clear symbols
- 24V batteries and a remote charging station so one battery can be charging while the other is in use

















1. Low Friction Casters

Low friction casters and dual precision bearings provide smooth, unsurpassed maneuverability that requires 60% LESS force to push/pull.

2. Integrated Base Cover

A unique integrated base cover provides excellent protection against wear and tear due to dirt, moisture or accidental damage.

3. Foldable Design

For added convenience, the I-Lift Series Lifts feature a unique foldable design for easy storage when not in use.

4. Redesigned Swivel Bar

The all new swivel bar design has double loop aluminum safety clips, a polyurethane wrapped protective cover and it can rotate 360° for maximum flexibility in positioning patients.

5. Premium Electronics

The I-Lift Power Plus Lift also includes a battery life indicator and digital service monitor. The digital monitoring light provides a convenient visual cue when service is due.

6. New Hand Pendant Design

I-Lift Series Lift models feature an ergonomic control which allows caregivers to operate the lift while at the patient's side with simple, easy to understand symbols.

Optional Digital I-Lift Scale

Easy to use with clearly labeled controls, an easily read display window and operating instructions printed on the face panel. Two-year limited warranty.

The Invacare® I-Lift™ Series Lifts offer Premium Features without the Premium Price.



Model # Height Base Casters Weight Weight Low Battery Emergency At Sling Hook Up Dimensions Capacity Alarm Controls ILIFTEM 22.2"-76.1" 4.6"H Clearance 118 lb. 450 lb. Audio Stop, Up 3.9" 37.4"W Open and Down 27.8"W Closed 46.6"L ILIFTEE 22.2"-76.1" 4.6"H Clearance 120 lb. 450 lb. Audio, Stop, Up 42.9"W Open and Down Visual 27.8"W Closed 46.6"L Charger Battery Charge Charger Limited Time Input Output Warranty ILIFTEM 24 V DC 100-240 29.5V DC 6 hrs. 4 Years Lift Rechargeable VAC 2 Years Electronics 2 w/Unit ILIFTEE 24 V DC 6 hrs. 100-240 29.5V DC 4 Years Lift

VAC

2 Years Electronics

Invacare Corporation

www.invacare.com

USA

One Invacare Way Elyria, Ohio 44035-4190 (800) 333-6900

Canada

© 2013

570 Matheson Blvd. E. Unit 8 Mississauga, Ontario L4Z 4G4 Canada (800) 668-5324

Invacare Corporation.
All rights reserved.
Trademarks are
identified by the
symbols ™ and ®.
All trademarks are
owned by or licensed

unless otherwise noted. Specifications are subject to change without notification. Form No. 13-222 130410

to Invacare Corporation



Complete instructions for lift usage and sling procedures are found in the Invacare Operating Manual. Only Invacare slings and accessories should be used on Invacare lifts.

Rechargeable

2 w/Unit