



Yes, you can.™

## Warranty

For all Lift Chairs

Chair frame	limited lifetime	limited lifetime	limited lifetime	limited lifetime
Fabric	1 year	1 year	1 year	1 year
Lift mechanism	limited lifetime	limited lifetime	limited lifetime	limited lifetime
Motor	3 years	3 years	3 years	3 years
Pendant	1 year	1 year	1 year	1 year
Body massage/heat	1 year	1 year	1 year	1 year
Labor	1 year	1 year	1 year	1 year

## Specifications

Two-Way Lift Chairs	LC2300	LC2700	LC2900
<b>Dimensions</b>			
Seat width	21.5"	21.5"	22.5"
Seat depth	21"	21"	21"
Floor-to-seat height	20"	20"	20"
Seat-to-back height	27"	27"	27.25"
Height of chair	43"	44.25"	44"
Overall depth	38"	37"	38"
Overall width	33.5"	35"	30.5"
Overall extension	60"	58.5"	60"
Recline angle	40°	37°	33°
<b>Weight</b>			
Weight of lift chair	144 lb.	145 lb.	136 lb.
Shipping weight	159 lb.	160 lb.	151 lb.

Three-Way Lift Chairs	LC3100	LC3300
<b>Dimensions</b>		
Seat width	21.5"	20"
Seat depth	21"	22"
Floor-to-seat height	19.5"	20"
Seat-to-back height	27.25"	25.5"
Height of chair	44"	42.5"
Overall depth	37.5"	37.5"
Overall width	33.5"	32"
Overall extension	67"	66"
Recline angle	60°	60°
<b>Weight</b>		
Weight of lift chair	131 lb.	133 lb.
Shipping weight	146 lb.	148 lb.

Three-Way Wall Hugger Lift Chairs	LC4100	LC4300	LC4300W	LC4500
<b>Dimensions</b>				
Seat width	21.5"	19.5"	25"	22.5"
Seat depth	20"	22"	21"	20"
Floor-to-seat height	19.5"	20.5"	21"	19.5"
Seat to-back height	26.5"	28.5"	27"	28"
Height of chair	41.5"	44"	44"	42"
Overall depth	37.5"	41"	39.5"	37.5"
Overall width	35"	34"	43"	34"
Overall extension	66.5"	70"	66"	66.5"
Recline angle	59°	55°	50°	60°
<b>Weight</b>				
Weight of lift chair	131 lb.	130 lb.	141 lb.	145 lb.
Shipping weight	146 lb.	145 lb.	156 lb.	160 lb.

Invacare Corporation  
www.invacare.com

### USA

One Invacare Way  
Elyria, Ohio  
44036-2125  
(800) 333-6900

### Canada

5970 Chedworth Way  
Mississauga, Ontario  
L5R 3T9 Canada  
(800) 668-5324

Invacare is a registered trademark and Yes, you can. is a trademark of Invacare Corporation. Dacron is a registered trademark of E. I. du Pont de Nemours and Company. Guardsman is a registered trademark of Lilly Industries, Inc. Wear-Dated is a registered trademark of Monsanto Company. Herculon is a registered trademark of Hercules Incorporated. Specifications are subject to change without notification.

Note: Contact Invacare Customer Service to determine availability of lift chair styles in your country of origin.

© 2000  
Invacare Corporation  
Form No. 00-247

Invacare® Lift Chairs



Yes, you can.™

Invacare's line of lift chairs is brand-new: re-engineered from top to bottom, from the inside out. And the result is outstanding. Designed and built by a leading manufacturer of motion furniture, these chairs feature a state-of-the-art lift mechanism and German-engineered motor to form the basis of quality chairs that will provide dependable performance, year after year. What's more, the Invacare Lift Chair line has expanded to include a wide chair and several space-saving Three-Way Wall Hugger models, with optional heat and massage. Finally, all Invacare Lift Chairs are attractively designed, with color and style choices that are sure to complement any consumer's décor for years to come.

## FEATURES

### Chair Upholstery

- A variety of styles and durable fabrics
- Fabrics meet California Bureau of Home Furnishings and Flammability I 17 standards and UFAC Class I
- Dacron® fiber in the arms, back and chaise helps retain shape after years of use
- Polyurethane foam in seat ensures comfort and support

### Chair Construction

- Engineered, multi-layered wood provides strength and long-term durability
- Frame is reinforced with mortise and tenon joints for strength
- Heavy-gauge, no-sag steel springs in seat and back provide durability and comfort
- Spring-loaded ottoman for safety

### Lift Mechanism and Actuator

- State-of-the-art lift mechanism is durable, reliable and dependable
- Powerful, low-voltage German-engineered motor is whisper-quiet, UL listed and CSA approved
- Welded steel base increases stability in elevated positions
- Key pivot points are engineered with bushings for dependability and smooth operation

### Freight and Handling

- Removable back allows for easy delivery and setup
- Compact carton requires less warehouse space and is stackable
- Free freight on any quantity

### General Features

- Weight capacity: 350 lb.
- Ergonomic hand pendants are easy to operate
- Options include:
  - Headrest covers and pillows in matching fabrics
  - Guardsman® fabric protection plan
  - Body Relaxation System: dual heat/massage (Three-Way Wall Hugger Lift Chairs only)
  - For an additional charge, any 2000 model fabric selection is available on any 2000 model, and any 3000 or 4000 model fabric is available on any 3000 or 4000 model

### Hand Pendants

- Pressure-sensitive controls
- Internationally recognized symbols
- 6-foot coil cord for convenience
- Three-Way Wall Hugger models are equipped with a safety lock system that prevents the chair from being raised or lowered when chair is unattended



Hand pendant for Two- and Three-Way Lift Chairs

Hand pendant for Three-Way Wall Hugger Lift Chairs

Hand pendant for Three-Way Wall Hugger Lift Chairs with Body Relaxation System



## Invacare® Two- and Three-Way Lift Chair Positions



Seated

Partially reclined



Fully reclined

Lift-out

## Two-Way Lift Chairs

Invacare Two-Way Lift Chairs are an excellent choice for the consumer who is looking for a high-quality chair at an excellent value. These chairs offer assistance and support by providing a partial reclining position for relaxation, as well as a lift-out position to aid in standing and exiting the chair.



## Invacare® LC2300

- T-back styling
- Quality comfort foam
- Three fabric selections

Contents

60% Polyester, 22% Acrylic, 15% Olefin, 3% Nylon



Beige  
1201-84



Copen Blue  
1201-97  
(shown)



Wine  
1201-77

### Invacare® LC2700

- Simulated headrest styling
- Plush fabric
- One magazine pocket
- High-density foam
- Four fabric selections

Content  
100% Nylon



Sandstone  
1046-84



Blue  
1046-95



Wine  
1046-78



Arbor Green  
1046-58  
(shown)



### Invacare® LC2900

- Contemporary coastal styling with natural oak armrests
- Pub back styling
- One magazine pocket
- High-density foam
- Five fabric selections

Content  
100% Olefin



Parchment  
Cream  
1992-82



Fawn Beige  
1992-83  
(shown)



Bark Variegated  
1992-84



Indigo Blue  
1992-97



Wine  
1992-78



## Three-Way Lift Chairs

Invacare Three-Way Lift Chairs are perfect for the consumer who is looking for a premium-quality chair that, in addition to the partial and lift-out positions, also offers a fully reclining position, making it easy to relax or even sleep.



### Invacare® LC3100

- Shirred back styling
- Two magazine pockets
- High-density foam
- Four fabric selections

Content  
100% Nylon



Sandstone  
1020-84  
(shown)



Blue  
1020-98



Burgundy  
1020-78



Arbor Green  
1020-58



### Invacare® LC3300

- Contemporary styling
- Pillow back
- Chaise pad
- Two magazine pockets
- High-density foam
- Four fabric selections

Contents  
79% Herculon®, 21% Polyester



Pecan  
1514-88  
(shown)



Sapphire  
1514-98



Ruby  
1514-78



Emerald  
1514-58

## Three-Way Wall Hugger Lift Chairs

Invacare's Three-Way Wall Hugger Lift Chairs are especially convenient for consumers who have limited space but would enjoy a chair that fully reclines. These chairs offer the same positions as the Three-Way Lift Chairs, but require only 4 to 6 inches of space between the back of the chair and the wall to fully recline, making them a perfect fit in any room.

### Body Relaxation System

Exclusively available on Invacare's Three-Way Wall Hugger Lift Chairs, the advanced Body Relaxation System offers soothing massage and comforting heat in the lumbar and seating area. An easy-to-use universal hand pendant controls the heat and massage feature along with the operation of the lift chair.

#### Invacare® LC4100

- Traditional button back
- Two magazine pockets
- High-density foam
- Safety lock system in hand pendant
- Optional Body Relaxation System: heat and massage
- Four fabric selections

Contents  
50% Polyester, 24% Acrylic, 24% Olefin, 2% Nylon



Taupe  
1741-84



Blue  
1741-95



Claret Mauve  
1741-77  
(shown)



Jade  
1741-57






LC4300W




### Invacare® LC4300/LC4300W

- Casual contemporary styling
- Pillow back
- Chaise pad
- Two magazine pockets
- High-density foam
- LC4300W available in 25" seat width
- Safety lock system in hand pendant
- Optional Body Relaxation System: heat and massage
- Four fabric selections  
Contents  
Wear-Dated® 58% Polyester, 26% Acrylic, 13% Cotton, 3% Nylon
- Four vinyl selections  
Contents  
98% Vinyl, 2% Polyurethane

Fabric

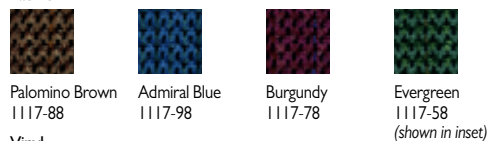

Vinyl

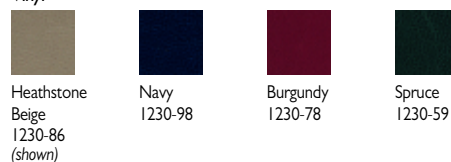

### Invacare® LC4500

- Traditional button back with skirt
- Tone-on-tone basketweave fabric pattern
- Pendant pocket
- High-density foam
- Safety lock system in hand pendant
- Optional Body Relaxation System: heat and massage
- Four fabric selections  
Contents  
52% Acrylic, 48% Polyester

