




Invacare® Premier Series  
Jasmine™ Full Body Lift

**Yes, you can.®**



Optimum functionality. The Jasmine lift has innovative features that provide safer resident handling within a 500 pound weight capacity. Caregivers can rely on safer operation and residents on maximum comfort during everyday transfers.

# Invacare® Premier Series Jasmine™ Full Body Lift

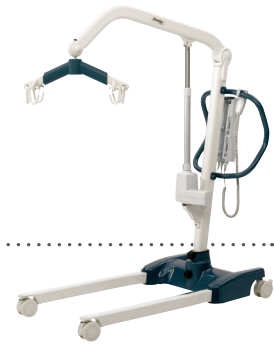

### Lifting Range

This lift expertly functions in more challenging situations with a rotating 360° 6-point hanger bar and 16.4" minimum height range, essential for assisting residents who may have fallen to the floor.



### Maneuverability

The push handle accommodates any staff member. And with a lighter-weight aluminum and steel construction, control in movement is straightforward and trouble-free.



### Ease For Caregivers

Jasmine lifts feature a battery-powered boom and wide-opening base, both controlled through the hand pendant



### Security

The stable spring-loaded hanger bar produces an audible tone when the sling is engaged. This assures the caregiver that the sling is properly attached.

## Jasmine™ Specifications

- Height at Sling Hook-up - MAX : 74.6 inches
- Height at Sling Hook-up - MIN : 16.4 inches
- Base Width OPEN : 46.7 inches
- Base Width CLOSED : 25.8 inches
- Base Height (Clearance) : 4.8 inches
- Base Length : 49 inches
- Caster Size (Front/Rear) : 4.0 inches
- Weight Capacity : 500 lbs/227 Kg.
- Weight Out of Carton : 98 lbs
- Battery : 24V DC (Rechargeable)
  - 2 Shipped with Unit
- Charger Input : 100-240V AC
- Charger Output/Charging Time : 29.5 DC Max 6 hours
- Audio Low Battery Alarm : Yes
- Motor Safety Devices : Anti-Entrapment
- Warranty Electronics : 3 Years On Lift, 2 Years on Electronics
- Emergency Stop Button : Yes
- Emergency Up & Down : Yes

## Jasmine™ Accessories

- Cradle
- XL Cradle



**Warning:** Before using the Patient Sling and Patient Lift, read and understand the Owner's Manual for proper operation and safety procedures.

- Scale : JLS5



**Invacare Continuing Care, Inc.**  
A Subsidiary of Invacare Corporation

[www.invacare-ccg.com](http://www.invacare-ccg.com)

[info\\_iccg@invacare.com](mailto:info_iccg@invacare.com)

### USA

One Invacare Way  
Elyria, Ohio 44035  
USA  
800-668-2337

### Canada

994 Hargrieve Road  
London, Ontario  
N6E 1P5 Canada  
800-668-2337

Trademarks are identified by the symbol™ and®.  
All trademarks are owned or licensed to Invacare Corporation or its subsidiaries unless otherwise noted.



**Yes, you can.®**