



Yes, you can:

Invacare® microAIR® MA80 & MA85 True Low Air Loss Mattresses




MA85
shown on CS7™ bed
with ThinkSoft® Positioning Device

These two (2) true low air loss mattress models help manage skin moisture to avoid maceration and maintain a dryer environment for healing pressure ulcers in the treatment of Stage I to Stage IV pressure ulcers.

Invacare® microAIR® MA80 & MA85 True Low Air Loss Mattresses



Body Interface Pressure Imaging Test

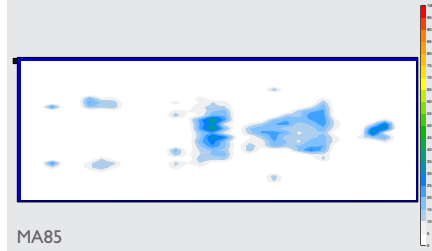


Test demonstrates contact measurement between key pressure points. Subject is 5'5" female, 110 lbs Simulation. Individual resident results will vary.



User-Friendly Control Unit

Lightweight unit with LED display and up down arrows for easy adjustment of comfort and therapy settings

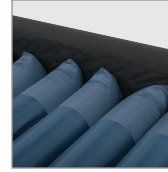

Quick-Connect Air Cells

Quick-connect air cells allow easy replacement of individual cells if necessary



90 degree Magnetic Hose Connector

MA85 features secure magnetic attachment simplifies set-up and provides quick CPR deflate



Wireless Auto Fowler

Automatically adds 25% more air to prevent bottoming out when mattress is articulated

Specifications

- Resident Risk Level : Medium-High
- Cover : Dual-layered, quilted, anti-microbial, fluid resistant vapor-permeable cover protects internal cells. Bottom cover is 1000 Denier nylon
- Limited Warranty :
1 year on mattress,
2 years on power unit
- Power Unit : 1275 LPM
- Electrical Power : 90V, 60Hz, 1A
- Internal cells : 70 Denier urethane coated nylon

Features

- Efficient real-time monitoring for correct pressure relief
- 2" safety foam base adds safety and comfort
- Static setting button
- Built-in carrying straps for easy transportation
- Settings retained during a power outage will resume once power is restored
- 1275 LPM blower inflates mattress in seconds
- MA85 includes 4 alternating pressure settings

	MA80	MA85	Bariatric			
			MA80B42	MA80B48	MA85B42	MA85B48
Width	36"	36"	42"	48"	42"	48"
Depth	10"	10"	10"	10"	10"	10"
Weight Capacity	350 lbs	350 lbs	1,000 lbs	1,000 lbs	1,000 lbs	1,000 lbs
Available with Raised Side Rails (RSR)	*	*	n/a	n/a	n/a	n/a

Invacare Continuing Care, Inc. www.invacare-ccg.com info_ccg@invacare.com

USA
One Invacare Way
Elyria, Ohio 44035
USA
800-668-2337

Canada
994 Hargrieve Road
London, Ontario
N6E 1P5 Canada
800-668-2337

For more information, contact your Invacare Continuing Care Sales Representative at 1.800.668.2337 or visit www.invacare-ccg.com.

© 2009 Invacare Continuing Care, Inc. All rights reserved. Trademarks are identified by the symbols ™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.

Form 08.66.9/09



Yes, you can.