



MKIV Joysticks

RII, TCHJOY, A, A+

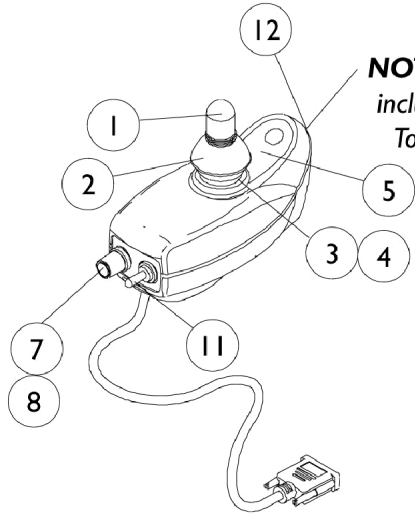


NOTE: RII 80 AMP. decal is not included with joystick. Used on Torque SP's after 10/31/01


1092610
RII Joystick

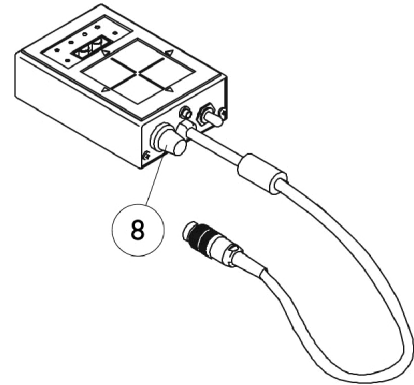

1095289
TCHJOY
Touchpad "A"
Joystick

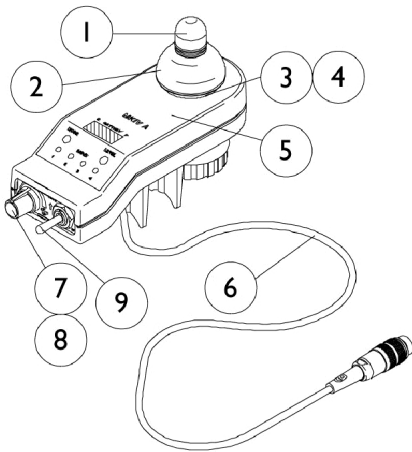

1065944
A Joystick

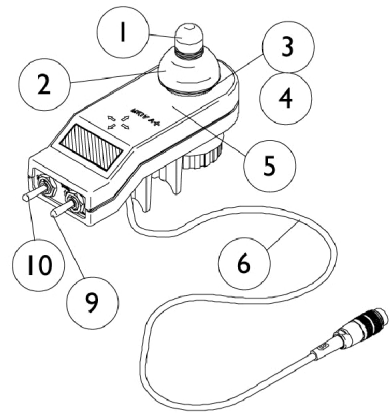

1038608
A+ Joystick

MKIV Joysticks RII, TCHJOY, A, A+

Item Number	Description	RII	TCHJOY	A	A+
See Pictorial	Joystick Assembly	1092610-NLA See Note(s) D	1095289-NLA	1065944-NLA	1038608-NLA
1	Knob, Joystick	1040217	N/A	1040217	1040217
2	Skirt, Joystick	1040218	N/A	1040218	1040218
3	Inductive, Joystick	1036764-NLA See Note(s) A	N/A	1036764-NLA See Note(s) B	1036764-NLA See Note(s) A
3	Inductive, Joystick	N/A	N/A	1104397-NLA See Note(s) C	N/A
4	Boot, Joystick	50000X364-NLA	N/A	50000X364-NLA	50000X364-NLA
5	Overlay, Joystick Decal	N/A	N/A	1040190-NLA	1036776-NLA
6	Cable Assembly	N/A	N/A	1036750-NLA See Note(s) B	1036750-NLA See Note(s) A
6	Cable Assembly	N/A	N/A	1104205-NLA See Note(s) C	N/A
7	Potentiometer, Speed	1040225-NLA See Note(s) A	N/A	1040225-NLA See Note(s) B	N/A
7	Potentiometer, Speed	N/A	N/A	1115518-NLA See Note(s) C	N/A
8	Knob, Potentiometer	1025516-NLA	1025516-NLA	1025516-NLA	N/A
9	Switch, On/Off Momentary	N/A	N/A	1040185-NLA See Note(s) B	1040185-NLA See Note(s) A
9	Switch, On/Off Momentary	N/A	N/A	1115523-NLA See Note(s) C	N/A
10	Switch, Program	N/A	N/A	N/A	1040186-NLA See Note(s) A
11	Switch, Toggle SPDT	1010042-NLA See Note(s) A	N/A	N/A	N/A
12	Decal, RII 80 AMP	1108692-NLA	N/A	N/A	N/A

NOTE: A - Item requires soldering
 B - Before 10/01/02 item requires soldering
 C - After 9/30/02 item has a non-soldering connector
 D - RII 80 amp. decal is not included with joystick. Used on Torque SP's after 10/31/01.