



Harnesses


As of 7/1/2009, switches are mercury free.

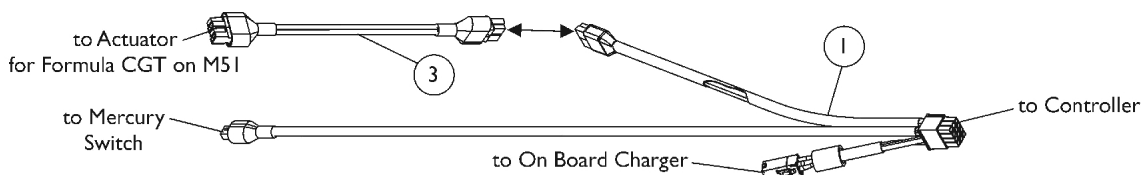


Tilt Actuator, Sensor, and Charger Harness (PTO Plus on M71/M91, M71jr Tilt chairs) & (Formula CGT on M51)

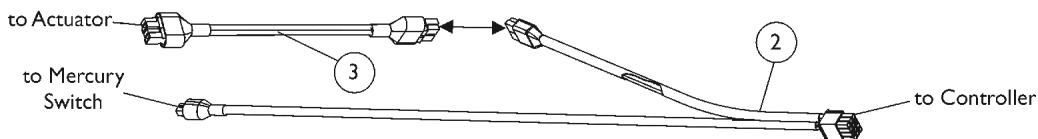


Tilt Actuator and Sensor Harness (Power Tiger w/Tilt), (Super Low Tilt & PTO Plus on TDX chairs), (Formula CG on TDX SP Before 6/19/07), (Formula CG on 3G Storm)

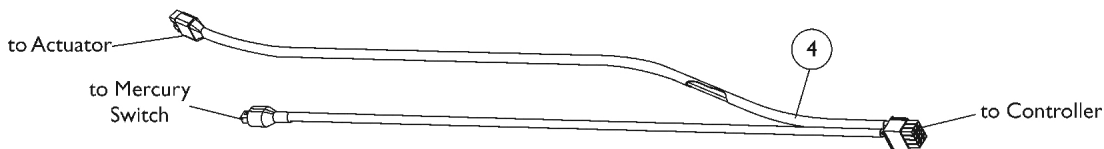


Tilt Actuator and Sensor Harness (Formula CG on TDX SP After 6/18/07) & (Formula CGT on TDX SC & TDX SI)

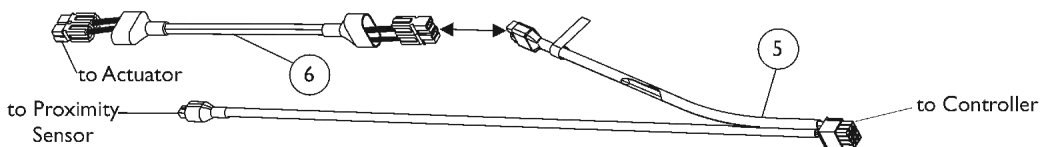

Elevate Actuator and Sensor Harness (Formula Elevate on TDX, Formula CGE on TDX SP, TDX SC & 3G Storm)

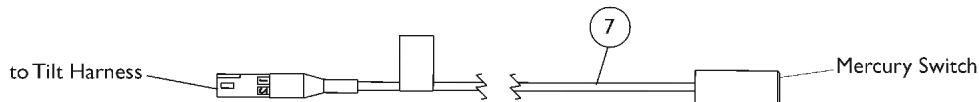

Drive Lockout Mercury Switch

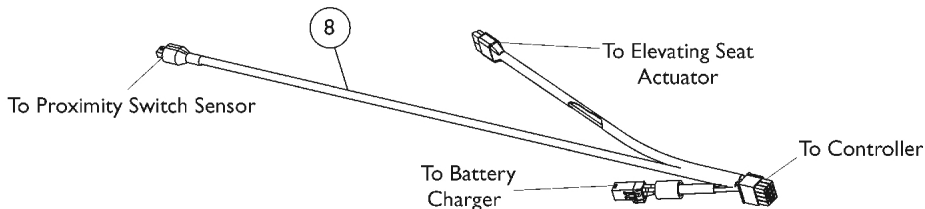

Elevate Actuator and Battery Charger Harness (Pronto M61)

Harnesses

As of 7/1/2009, switches are mercury free.

Item Number	Note	Part Number	Description	Quantity Per	
				Package	Assembly
1	A	1140124	Harness, Tilt Actuator with Sensor and Charger	1	1
2	B	1140101	Harness, Tilt Actuator and Sensor	1	1
3	C	1145414	Cable, Actuator Jumper Adapter (MK6i)	1	1
3	G	1101400	Harness, ELRPW Power Legrest - Right	1	1
4	D,D1	1150961	Harness, Tilt Actuator and Sensor	1	1
5	E,F	1140100	Harness, Elevate Actuator and Sensor	1	1
6	F	1142207	Cable, Actuator Jumper Adapter R2T (MK6i)	1	1
7	H	1140099	Switch, Drive Lockout Mercury Free Ball (MK6i) (50" L)	1	1
8	I	1127260	Harness, Actuator Sensor and Battery Charger	1	1

NOTE: A - Chairs with (PTO Plus on Pronto M71/M91), (M71jr tilt) and (Formula CGT on M51) models.
 B - Chairs with (Manual Tilt on Power Tiger), (Super Low Tilt & PTO Plus on TDx), (Formula CG on TDx SP before 6/19/07) and (Formula CG w/ a Conventional Tilt Actuator Motor on 3G Storm and 4-Pole Motors).
 C - Used for TDx SP models made before 6/19/07. If the customer wants to eliminate the need for the adapter cable, order item #4 "Harness, Tilt Actuator and Sensor" dated after 6/18/07
 D - After 6/18/07 a new Tilt Actuator and Sensor Harness was created to eliminate the need for item #3 "Cable, Actuator Jumper Adapter". The Harness is backward compatible for chairs made before 6/19/07 if the customer wants to eliminate the adapter cable.
 D1 - Used for a Conventional Tilt Actuator Motor only on TDx SP-CG (Tilt Only) models equipped with 4-Pole Motors. Used for Formula CGT on TDx SC and TDx SI models also.
 E - Chairs with (Formula Elevate on TDx), (Formula CGE on TDx SP, TDx SC and 3G Storm).
 F - The Actuator Jumper Adapter R2T Cable is only used in connection with item #5 Elevate Actuator and Sensor Harness for Formula Elevate on TDx
 G - Part # 1101400 traditionally used for ELRPW Legrest is also used for the Formula CG Tilt Actuator Motor on Pronto M51 bases. This Harness is also used for 3G Storm Series with 4-Pole Motors and a Conventional Tilt Actuator Motor.
 H - Mercury Switch is used for Conventional Actuator Motors only.
 I - Elevate Actuator and Battery Charger Harness (Pronto M61).