


Compressor and Mounting Hardware Platinum 5 Models

NOTE: Compressors are interchangeable. However, Capacitor and Bracket Restraint must also be replaced at the same time.

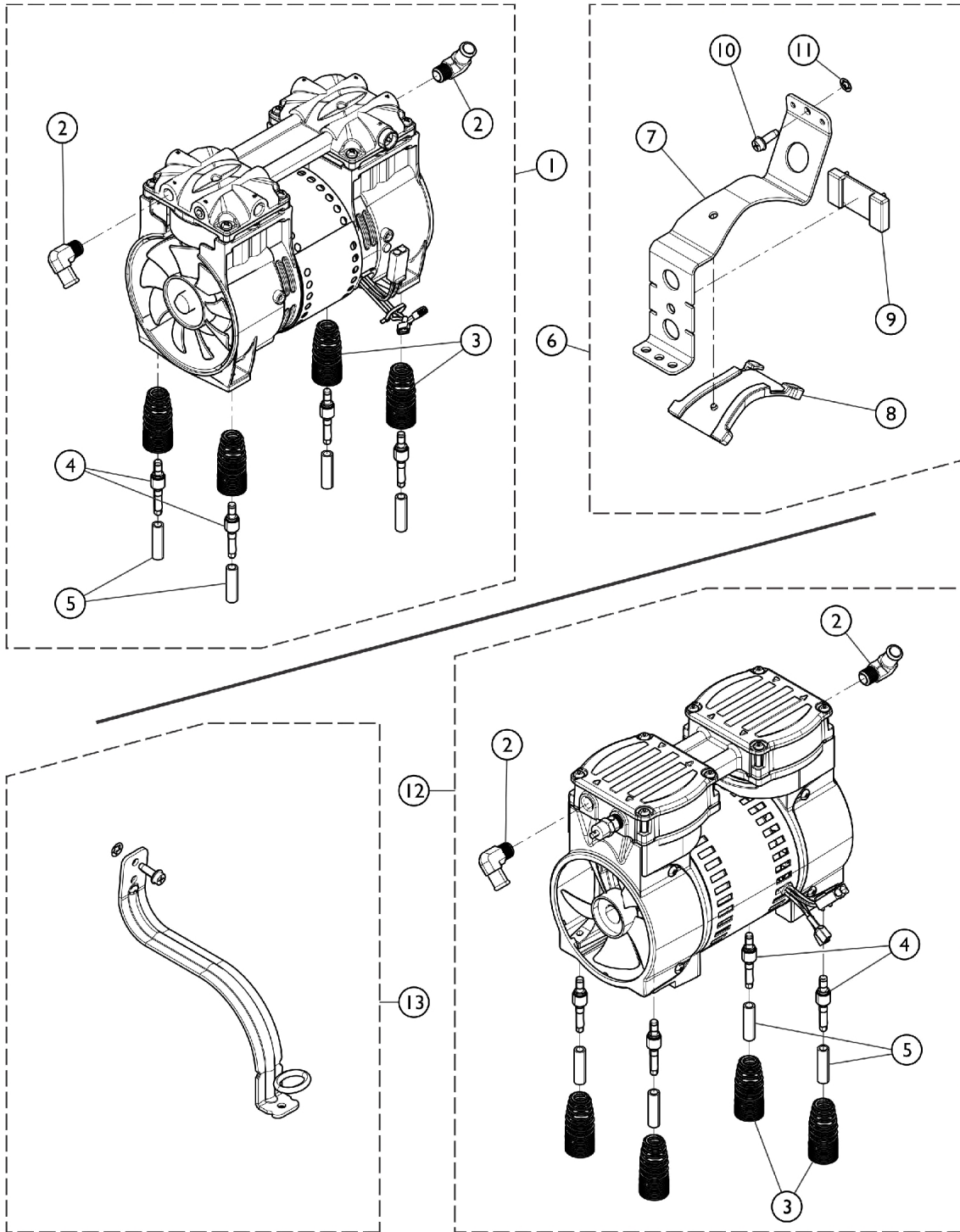

Compressor and Mounting Hardware

Platinum 5 Models

NOTE: Compressors are interchangeable. However, Capacitor and Bracket Restraint must also be replaced at the same time.

Item Number	Note	Part Number	Description	Quantity Per				
				Package	Assembly			
1	A	1157138	Kit, Compressor Complete Assembly - Thomas	1	1			
			Compressor with Mounting Hardware - Thomas	1	1			
			NLA - Tech Only - Capacitor, 15 MFD	1	1			
			Tie Wrap (8" L)	2	1			
		1186230	Kit, Compressor Restraint Bracket w/ Hardware	1	1			
			Bracket, Restraint, Compressor	1	1			
			NLA - Bumper, Bracket	1	1			
			Bumper, Restraint Bracket, Compressor	1	1			
			NLA - Screw, Phillips Hex Washer Head w/ Patch (1/4-20 x 3/4")	1	1			
			Retainer, Push (1/4")	1	1			
1	B	1106690E-NLA	NLA - WARR RPLMT - Compressor Assembly - Thomas (5LX Plat.)	1	1			
			2	1108819	Fitting, Elbow (1/2" Beaded Hose x 1/4" NPT)	1	2	
			3	1105504-NLA	NLA - Spring, Compressor Mount	1	4	
			3	1150581	Spring, Compressor Mount	1	4	
			4	1106446-NLA	NLA - Pin, Spring Retainer	1	4	
			4	1150580	Pin, Spring Retainer	1	4	
			5	1107411	Glide, Pin	1	4	
			6	1186230	Kit, Compressor Restraint Bracket w/ Hardware	1	1	
			7		Bracket, Restraint, Compressor	1	1	
			8		NLA - Bumper, Bracket	1	1	
			9		Bumper, Restraint Bracket, Compressor	1	1	
10		NLA - Screw, Phillips Hex Washer Head w/ Patch (1/4-20 x 3/4")	1	1				
11	C,D,B	1154510	Retainer, Push (1/4")	1	1			
			Kit, Compressor Rebuild w/ Trigon Seals & Instruction Sheet - IRC5LX/IRC5LXO2 (Thomas)	1	1			
			Compressor Rebuild Kit w/ Trigon Seals - IRC5LX/IRC5LXO2 (Thomas)	1	1			
			Instruction Sheet, Thomas Rebuild 2660-2650	1	1			
			D,E	1068698-NLA	NLA - Tech Only - Capacitor, 15 MFD	1	1	
			D	2001682	Tie Wrap (8" L)	1	2	
			12	F	1125371-NLA	NLA - Compressor with Mounting Hardware - Gast	1	1
			12	G	1125371E-NLA	NLA - Compressor Assembly, Warranty Replacement - Gast (5LX Plat.)	1	1
			D,H	1124028	Kit, Rebuild for Gast Compressors	1	1	
					Kit, Rebuild	1	1	
					Instruction Sheet - Compressor Rebuild	1	1	
D,I,J	1131250KIT	Kit, Capacitor (10 MFD)	1	1				
		NLA - Tech Only - Capacitor, 10 MFD	1	1				
		Tie Wrap (17-1/2" L)	1	1				
		Tie Wrap (8" L)	1	2				
		Instruction Sheet, Capacitor Replacement	1	1				
13	J	1121244	Kit, Compressor Restraint Bracket (Gast)	1	1			

Compressor and Mounting Hardware

Platinum 5 Models

NOTE: Compressors are interchangeable. However, Capacitor and Bracket Restraint must also be replaced at the same time.

NOTE: A - Kit includes items 2 - 6 and 1186230. Must use with Thomas Compressors.
B - When replacing or servicing a Compressor, the Capacitor (1068698) and two Tie Wraps (2001682) must also be replaced.
C - Kit includes Retainer Plates and Trigon Seals to enable rebuilding of both Trigon and pre-Trigon configurations.
D - Not shown.
E - Must use with Thomas Compressors.
F - Item is No Longer Available, refer to item 1 (1157138) as its replacement.
G - When replacing or servicing a Compressor, the Capacitor (1131250) must also be replaced.
H - Includes Instruction Sheet, (2) Head O-Rings, (4) Leaf Valves, (4) Self Tapping Screws, (4) Retainer Screws, and (2) Piston Cups.
I - Includes Capacitor, Tie Wraps, and Installation Instructions.
J - Must use with Gast Compressors.