



Yes, you can.®

Pronto M41™



Anything else is a compromise.™



Yes, you can.®

The Invacare® Pronto® M41™ is the latest offering in the sleek and versatile Pronto® series of Power wheelchairs. The M41™ features the latest technology in electronics and in-line motors. This combination provides increased efficiency, battery range and an increased top speed. The M41™ also features a center wheel drive platform, simplified SureStep®, and six wheels on the ground providing stability and maneuverability similar to other Pronto® models. All of these features are packed into a simple, clean-looking, and easy-to-service compact frame that provides mobility with style for people on the go!




Pronto® M41™ with Fold-Down Seat (Standard)



Pronto® M41™ with Semi-Recline Van Seat (optional)

Colors Available:

Candy Red

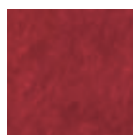

Deep Blue



Pronto M41 Features:

- Compact size for a tight turning radius
- High efficiency in-line motors for greater speed and longer battery range
- SureStep® suspension with six wheels on the ground for smooth driving and optimum stability
- Charger connection conveniently located on joystick
- Stylish design
- 3-amp off-board charger
- Comfortable and lightweight fold down seat with width adjustable and fixed height arms (standard)
- Semi-recline van seat with width and height adjustable arms (optional)
- Auto style positioning strap
- U1 Batteries
- 300 lb. weight capacity

Pronto® M41™ Specifications

Overall Dimensions (measured with fold down seat)		
Length	34"	
Width (including joystick)	24"	
Height	33"	
Seat-to-Floor	19"	
Seat Dimensions (fold down seat) (semi-recline seat)		
Width	18"	18"
Depth	17"	16" - 18"
Seat Angle	92°	80° - 115°
Weight capacity	300 lb.	
Wheel Sizes		
Front/Rear Casters	6" x 2"	
Drive Wheels (Flat Free)	10.5" x 3.5"	
Maximum Speed	5 mph	
Motor	250 watts	
Incline Capability	7° grade	
Turning Radius		
Front	17.75"	
Rear	19.5"	
Weight		
with U1 Batteries	161 lb.	
without batteries and foot plate	110 lb.	
Ground Clearance	2"	
Battery	U1	
Battery Charger	3 amp off-board	
Per-Charge Battery Range	Up to 17 miles	
Limited warranties		
Base frame	36 months	
Seat frame	36 months	
Electronics	12 months	
Motor/Gearbox	12 months	

Invacare Corporation
www.invacare.com

USA
One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada
570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

Invacare, the Blue Medallion Design, Yes, you can., Pronto and SureStep are registered trademarks of Invacare Corporation. M41 and Anything else is a compromise, are Trademarks of Invacare Corporation. Specifications are subject to change without notification.

© 2006
Invacare Corporation
Form No. 06-033



"Invacare says 'Yes, you can.' and I couldn't have said it better myself."